



# Composi-Tight® 3D Fusion™

## NEW! Full Curve Non-Stick Matrices

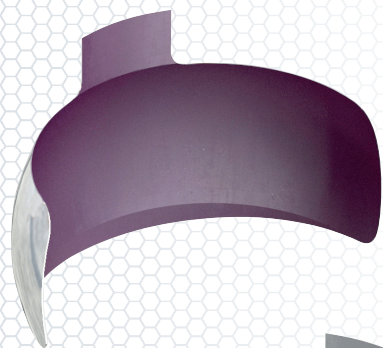

**Wrap it up** – Full Curve bands wrap further around the tooth, out of the way, making ring placement easier. Still dead-soft and a skinny 0.0016" thin, they're the perfect choice.


**Keep tabs on 'em** – An occlusal placement Grab-Tab™ gives you excellent control with any placement instrument you choose. Just bend it over the adjacent tooth after you've positioned the band.

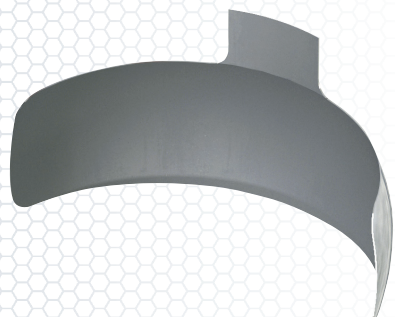



**Go deep or wide** – Improved subgingival extensions and added wrap-around length mean improved results on those more challenging restorations.



**Easy in** – Even though 3D Fusion Full Curve Matrices are made from the same ultra-thin dead-soft stainless steel as our other non-stick bands, their advanced curved form improves structural rigidity, making them our easiest-to-place band ever.

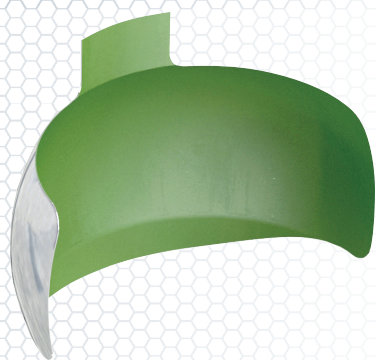

**Easy out** – 3D Fusion Full Curve matrices wouldn't be complete without the Slick Bands coating invented by Garrison. It reduces the adhesion of bonding agents and composites by a massive 92%. Regardless of your materials, these bands will slide right out!

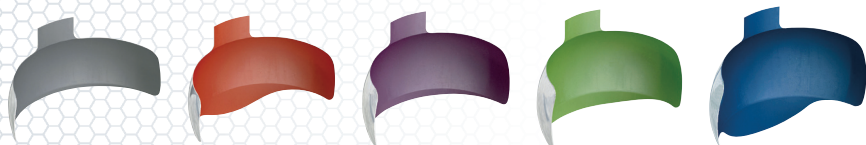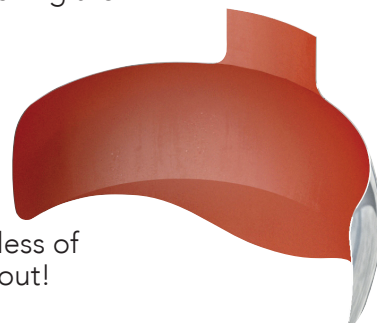

### 3D Fusion™ Bands

#### Description

Fusion Band Kit: 5 sizes  
Fusion Band Kit: 5 sizes

#### Item

FXB01  
FXB02

#### Qty.

420  
210

#### Refills

Premolar 4.4 mm  
Premolar w/extension 6.0 mm  
Molar 5.6 mm  
Tall Molar 6.6 mm  
Tall Molar w/extension 8.7 mm

#### Item

FX100-M  
FX150-M  
FX175-M  
FX200-M  
FX300-M

#### Standard

#### Qty.

50  
30  
50  
50  
30

#### Bulk

#### Item

FX100  
FX150  
FX175  
FX200  
FX300

#### Qty.

100  
60  
100  
100  
60



**Garrison**  
Dental Solutions

# Composi-Tight® 3D Fusion™

## Ultra Adaptive Wedges

### Combining Soft-Face™ adaptive materials,

a firm inner core and advanced mechanical features to produce a wedge that truly works. Every time.

**Get a grip** The universal gripping block gives you absolute control of the wedging process whether you're using a cotton pliers, hemostat or "pin" tweezers. Use what you like – 3D Fusion has you covered.



**A real softie** 3D Fusion's Soft-Face over-mold allows the wedge to do what no other wedge can truly do – actually adapt to interproximal irregularities.

#### Tough on the inside

The firm inner skeleton allows for easy insertion and tooth separation just like a traditional wedge.

#### No pain – all gain


The upturned ski-like tip and soft radius edges glide over the rubber dam or papilla for trauma free insertion.

#### Keep it clean (and dry!)

3D Fusion's patent pending design easily follows the contours of the teeth sealing things up nice and tight while preventing overhangs.

#### Fins with serious function

Soft retentive fins smoothly fold down during wedge insertion and then spring back when clear of the interproximal space. Wedge back-out is a thing of the past.

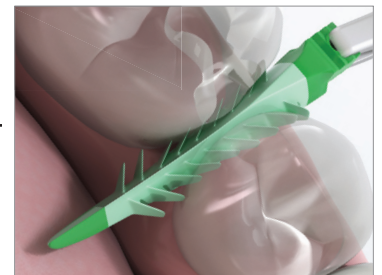

### 3D Fusion™ Wedge

Description	Item	Qty.
3D Fusion Wedge Kit	FXK4	400
3D Fusion Wedge Kit	FXK4-M	200

Refills	Item	Standard	Item	Bulk
		Qty.		Qty.
X-Small Yellow	FXYL-M	50	FXYL	100
Small Blue	FXBL-M	50	FXBL	100
Medium Orange	FXOR-M	50	FXOR	100
Large Green	FXGR-M	50	FXGR	100



For more information go to [www.garrisdental.com](http://www.garrisdental.com). Be sure to check out the website to see the full line of innovative Garrison products and our 6-month money-back guarantee.

**Garrison**  
Dental Solutions

150 DeWitt Lane  
Spring Lake, MI USA  
49456

To order: Call toll free **888.437.0032**  
or online at [www.garrisdental.com](http://www.garrisdental.com)



Mention Code  
**INROW-0816FXB**